

Autoclave Vertical



Double wall construction with inner S.S. 304. all these sheets are argon welded & outer chamber made of high quality S.S. 304/M. S sheet duly powder coated. Electrically heated by means of I.S.I. quality immersion heaters with working pressure 15 PSI. Hydraulically dye pressed Lid from S.S. Plate with. easy pedal lifting device. having pressure gauge. two Nos. safety valve. quick release coupling and handle with heater cover stand. Water level arrangement cut off buzzer system and S.S. perforated basket supplied along with unit on free of cost.

Inner S.S & Outer M.S


Model No.	Dia x Height	Rating	Volume
EAUC-A1	250 x 450mm	1.5KW	22Ltr
EAUC-A2	300 x 500mm	2.0KW	35Ltr
EAUC-A3	350 x 550mm	3.0KW	53Ltr
EAUC-A4	450 x 600mm	4.0KW	90Ltr
EAUC-A5	550 x 750mm	6.0KW	175Ltr

Inner & Outer S.S

EAUC-S1	250 x 450mm	1.5KW	22Ltr
EAUC-S2	300 x 500mm	2.0KW	35Ltr
EAUC-S3	350 x 550mm	3.0KW	53Ltr
EAUC-S4	450 x 600mm	4.0KW	90Ltr
EAUC-S5	550 x 750mm	6.0KW	175Ltr

OPTIONAL

- a) Digital Temperature controller cum indicator (instead of Energy Regulator)
- b) Digital Timer
- c) Auto Pressure cut off system
- d) Auto Vacuum Breaker
- e) Auto Purging & Exhaust (Solenoid Valve)
- f) Water level glass indicator

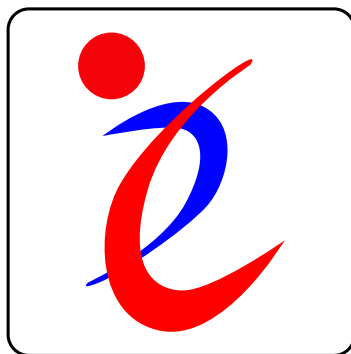

ELS Instruments

**Call Us****+91 94488 27690****Reach Us at****enquiry@elsinstruments.com**